

# Electric Breast Pump

HI-BEBE SUPER





### Portability

#### Rechargeable Battery

Don't need to bring power adaptor for expression up to 150 min. Feel the freedom of electric cable.

#### Light Weight and Small Size

Babies are tiny but their necessities are not. Small size and light weight help mothers carry the pumps anywhere with baby care product.

### Performance

#### Various Levels for Expression & Massage Mode

Find your own style using various expression and massage levels

(Expression : 16 steps, 6 cycles / Massage : 16 steps, 3 cycles)

BT-150 is designed desktop type in detail to control accurately and it can be used for hospital purpose. Save your time from powerful & accurate expression mode.

### Convenience

#### Nursing Night Lamp

Don't need another light if you turn on 'Nursing Night Lamp' and it creates a warming atmosphere.

#### Program

Make your own comfortable pattern using program function.

You can make 8 Program IDs & 8 Program sequences each IDs.

#### Mobile Application

- Manage and control using mobile application
- Speech recognition
- Compatible with Android & iPhone
- Setting program
- Control expression, massage mode and nursing night lamp level

### Hygiene

#### Closed Pumping System

Generated air from the pump blocked by diaphragm preventing infection. Also, this structure prevents backflow of milk into the tube and the pump.

#### BPA Free / Latex Free

All the parts contact the milk are made of BPA And Latex free material.



# BISTOS Technical Specification

Model	BT-150 (Desktop Type)	
Product Pictures		
Order code	150S	150D
Funnel & Bottle set	1 set	2 set
Vacuum Level	-20 to -290 mmHg	
Vacuum Control	16 Step	
Expression Cycle	35 to 60 Cycle / Min (Expression)	
Massage Cycle	70 to 90 Cycle / Min (Massage)	
Cycle Control	6 Step (Expression) / 3 Step (Massage)	
Battery Type (Capacity)	7.4 V Li-ion (2200 mA)	
Battery Run Time	120 Min	
Battery Charging Time	150 Min	
Funnel Size	24 mm (Option: 21 / 27 / 32 mm)	
Display	60 mm diameter Mono LCD	
Mode	Expression, Massage, Memory / Program	
Optional Accessories	Carrying Bag, Car Charging Cable, Funnel & Bottle set for double expression * Nursing Night Lamp, Wireless Mobile App	
Dimension (H x W x D)	190 x 190 x 79 mm	
Weight	1,100 g	
Pumping System	Closed System (Backflow Prevention)	
Auto Shut off	30 Min	
Materials	BPA Free / Latex Free	
Warranty	1 year (Main unit)	
Certificates	FDA, CE, ISO13485, GMP	
Origin	Rep. of Korea	




Carrying Bag



Nursing Night Lamp



Funnel Size 21 / 24 / 27 / 32 mm

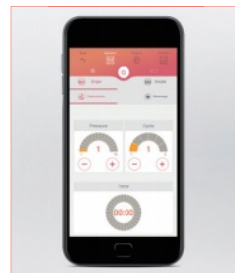

Mobile Application

## BIO SIGNAL TOTAL SOLUTION

### Bistos Co., Ltd. (Headquarter)

7th Fl., A Bldg., Woolim Lions Valley 5-cha, 302,  
Galmachi-ro, Jungwon-gu, Seongnam-si,  
Gyeonggi-do, Korea  
Tel. +82 (31) 750 0340  
Fax. +82 (31) 750 0344

bistos@bistos.co.kr  
www.bistos.co.kr



### Bistos Co., Ltd.

Bistos America Inc. (US Branch)  
22941, Triton Way, Suite #242, Laguna Hills, CA., USA 92653  
Tel. +1-949 614 8745 Fax. +1-949 614 8745  
Site. www.bistosamerica.com  
E-mail. info@bistosamericainc.com  
bistos@bistos.co.kr  
www.bistos.co.kr